


Compact laboratory balance KERN EWJ



**Promocija**  
**KEFO**  
Do kraja  
2019.

**BA**  
QUALITY

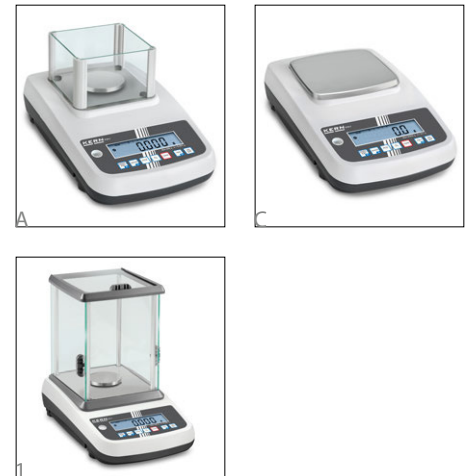

High-quality precision balance with automatic internal adjustment, also with EC type approval [M]

**Features**

- Easy to use: All primary functions have their own key on the keypad
- Automatic internal adjustment , time-controlled every 2 h, guarantees high degree of accuracy and makes the balance independent of its location
- Capacity display: A bar lights up to show how much of the weighing capacity is still available
- KERN EWJ/EWJ-H/EWJ-M: USB interface for transferring weighing data to the PC, printer etc., possibly in conjunction with KERN DBS-A02
- **Draft shield** standard on models KERN EWJ 300-3, EWJ 600-2SM, EWJ 600-2M, weighing space WxDxH 160x145x80 mm
- 1 KERN EWJ 300-3H: Large glass draught shield with 3 sliding doors for easy access to the items being weighed. Weighing space WxDxH 175x155x217 mm

- Protective working cover included with delivery

**Technical data**

- Large backlit LCD display, digit height 16,5 mm
- Dimensions weighing surface  
A 80 mm  
B 120 mm, see larger picture  
C WxD 155x145 mm
- Overall dimensions WxDxH  
220x315x90 mm (without draught shield)  
220x340x321 mm (incl. draught shield)
- Permissible ambient temperature  
KERN EWJ: 15 °C/35 °C  
KERN EWJ-M: 15 °C/30 °C


**Accessories**

- Protective working cover , scope of delivery: 5 items, KERN EWJ-A04S05
- Rechargeable battery pack internal , operating time up to 15 h without backlight, charging time approx. 4 h, KERN KFB-A01
- KERN EWJ-/EWJ-H/EWJ-M: USB Accessory Kit for bi-directional data exchange between balance/moisture analyser and computer. Scope of delivery: USB cable, driver, Software BalanceConnection, KERN DBS-A02
- RS-232/WiFi adapter for wireless connection to networks and WiFi capable devices, such as tablets, laptops or smartphones, KERN YKI-03
- RS-232/Ethernet adapter for connection to an IP-based Ethernet network, KERN YKI-01
- Further details, plenty of further accessories and suitable printers see Accessories

**STANDARD**



**OPTION**

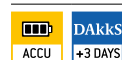

**FACTORY**



Model	Weighing capacity [Max] g	Readability [d] g	Verification value [e] g	Minimal load [Min] g	Linearity g	Weighing plate	Price	Option			
								Verification	DAkks Calibr. Certificate		
KERN								M KERN		DAkks KERN	
EWJ 300-3	300	0,001	-	-	± 0,005	A	€327	-		963-127	
EWJ 300-3H	300	0,001	-	-	± 0,005	A	€327	-		963-127	
EWJ 3000-2	3000	0,01	-	-	± 0,05	B	€312	-		963-127	
Note: For applications that require verification, please order verification on at the same time, initial verification at a later date is not possible. Verification at the factory, we need to know the full address of the location of use.											
EWJ 600-2SM	600	0,01	0,1	0,5	± 0,03	B	€265	965-216		963-127	
EWJ 600-2M	600	0,01	0,1	0,5	± 0,03	B	€300	965-216		963-127	
EWJ 6000-1SM	6000	0,1	1	5	± 0,3	C	€257	965-217		963-128	
EWJ 6000-1M	6000	0,1	1	5	± 0,3	C	€265	965-217		963-128	

## Pictograms

<b>CAL INT</b>	<b>Internal adjusting:</b> Quick setting up of the balance's accuracy with internal adjusting weight (motor driven)	<b>KCP PROTOCOL</b>	<b>KERN Communication Protocol (KCP):</b> It is a standardized interface command set for KERN balances and other instruments, which allows retrieving and controlling all relevant parameters and functions of the device. KERN devices featuring KCP are thus easily integrated with computers, industrial controllers and other digital systems	<b>IP</b>	<b>Protection against dust and water splashes IPx:</b> The type of protection is shown in the pictogram.
<b>CAL EXT</b>	<b>Adjusting program CAL:</b> For quick setting up of the balance's accuracy. External adjusting weight required	<b>GLP INTERN</b>	<b>GLP/ISO log:</b> The balance displays serial number, user ID, weight, date and time, regardless of a printer connection	<b>INOX</b>	<b>Stainless steel:</b> The balance is protected against corrosion
<b>EASY-T</b>	<b>Easy Touch:</b> Suitable for the connection, data transmission and control through PC, tablet or smartphone	<b>GLP PRINTER</b>	<b>GLP/ISO log:</b> With weight, date and time. Only with KERN printers	<b>UNDER</b>	<b>Suspended weighing:</b> Load support with hook on the underside of the balance
<b>MEMORY</b>	<b>Memory:</b> Balance memory capacity, e.g. for article data, weighing data, tare weights, PLU etc.	<b>PCS</b>	<b>Piece counting:</b> Reference quantities selectable. Display can be switched from piece to weight	<b>BATT</b>	<b>Battery operation:</b> Ready for battery operation. The battery type is specified for each device
<b>ALIBI</b>	<b>Alibi memory:</b> Secure, electronic archiving of weighing results, complying with the 2014/31/EU standard.	<b>GLP PRINTER</b>	<b>GLP/ISO log:</b> With weight, date and time. Only with KERN printers	<b>ACCU</b>	<b>Rechargeable battery pack:</b> Rechargeable set
<b>RS 232</b>	<b>Data interface RS-232:</b> To connect the balance to a printer, PC or network	<b>RECIPE A</b>	<b>Recipe level A:</b> The weights of the recipe ingredients can be added together and the total weight of the recipe can be printed out	<b>MULTI</b>	<b>Universal mains adapter:</b> with universal input and optional input socket adapters for A) EU, CH; B) EU, CH, GB, USA; C) EU, CH, GB, USA, AUS
<b>RS 485</b>	<b>RS-485 data interface:</b> To connect the balance to a printer, PC or other peripherals. Suitable for data transfer over large distances. Network in bus topology is possible	<b>RECIPE B</b>	<b>Recipe level B:</b> Internal memory for complete recipes with name and target value of the recipe ingredients. User guidance through display	<b>230 V</b>	<b>Mains adapter:</b> 230V/50Hz in standard version for EU. On request GB, USA or AUS version available
<b>USB</b>	<b>USB data interface:</b> To connect the balance to a printer, PC or other peripherals	<b>RECIPE C</b>	<b>Recipe level C:</b> Internal memory for complete recipes with name and target value of the recipe ingredients. User guidance through display, multiplier function, adjustment of recipe when dosages are exceeded or barcode recognition	<b>230 V</b>	<b>Power supply:</b> Integrated in balance. 230V/50Hz standard EU. More standards e.g. GB, USA or AUS on request
<b>BT</b>	<b>Bluetooth* data interface:</b> To transfer data from the balance to a printer, PC or other peripherals	<b>SUM A</b>	<b>Totalising level A:</b> The weights of similar items can be added together and the total can be printed out	<b>DMS</b>	<b>Weighing principle: Strain gauges</b> Electrical resistor on an elastic deforming body
<b>WLAN</b>	<b>WLAN data interface:</b> To transfer data from the balance to a printer, PC or other peripherals	<b>PERCENT</b>	<b>Percentage determination:</b> Determining the deviation in % from the target value (100 %)	<b>T-FORK</b>	<b>Weighing principle: Tuning fork:</b> A resonating body is electromagnetically excited, causing it to oscillate
<b>SWITCH</b>	<b>Control outputs (optocoupler, digital I/O):</b> To connect relays, signal lamps, valves, etc.	<b>UNIT</b>	<b>Weighing units:</b> Can be switched to e.g. nonmetric units at the touch of a key. See balance model. Please refer to KERN's website for more details	<b>FORCE</b>	<b>Weighing principle: Electromagnetic force compensation</b> Coil inside a permanent magnet. For the most accurate weighings
<b>ANALOG</b>	<b>Analogue interface:</b> to connect a suitable peripheral device for analogue processing of the measurements	<b>TOL</b>	<b>Weighing with tolerance range:</b> (Checkweighing) Upper and lower limiting can be programmed individually, e.g. for sorting and dosing. The process is supported by an audible or visual signal, see the relevant model	<b>SC TECH</b>	<b>Weighing principle: Single cell technology:</b> Advanced version of the force compensation principle with the highest level of precision
<b>DUAL</b>	<b>Interface for second balance:</b> For direct connection of a second balance	<b>MOVE</b>	<b>Hold function:</b> (Animal weighing program) When the weighing conditions are unstable, a stable weight is calculated as an average value	<b>+3 DAYS</b>	<b>Verification possible:</b> The time required for verification is specified in the pictogram
<b>LAN</b>	<b>Network interface:</b> For connecting the scale to an Ethernet network			<b>DAKkS +3 DAYS</b>	<b>DAKkS calibration possible:</b> The time required for DAKkS calibration is shown in days in the pictogram
<b>RC</b>	<b>Wireless data transfer:</b> between the weighing unit and the evaluation unit using an integrated radio module			<b>1 DAY</b>	<b>Package shipment:</b> The time required for internal shipping preparations is shown in days in the pictogram

\*The Bluetooth® word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. and any use of such marks by KERN & SOHN GmbH is under license. Other trademarks and trade names are those of their respective owners.

## KERN – Precision is our business

To ensure the high precision of your balance KERN offers you the the appropriate test weight in the international OIML error limit classes E1-M3 from 1 mg - 2500 kg. In combination with a DAKkS calibration certificate the best pre-requisite for proper balance calibration.

The KERN DAKkS calibration laboratory today is one of the most modern and best-equipped DAKkS calibration laboratories for balances, test weights and force-measurement in Europe.

Thanks to the high level of automation, we can carry out DAKkS calibration of balances, test weights and force-measuring devices 24 hours a day, 7 days a week.

### Range of services:

- DAKkS calibration of balances with a maximum load of up to 50 t
- DAKkS calibration of weights in the range of 1 mg – 2500 kg
- Volume determination and measuring of magnetic susceptibility (magnetic characteristics) for test weights
- Database supported management of checking equipment and reminder service
- Calibration of force-measuring devices
- DAKkS calibration certificates in the following languages DE, GB, FR, IT, ES, NL, PL
- Conformity evaluation and reverification of balances and test weights

## Your KERN specialist dealer:

